




LABORATORY WORKTOPS
INTERIOR PANELS



Think Trespa

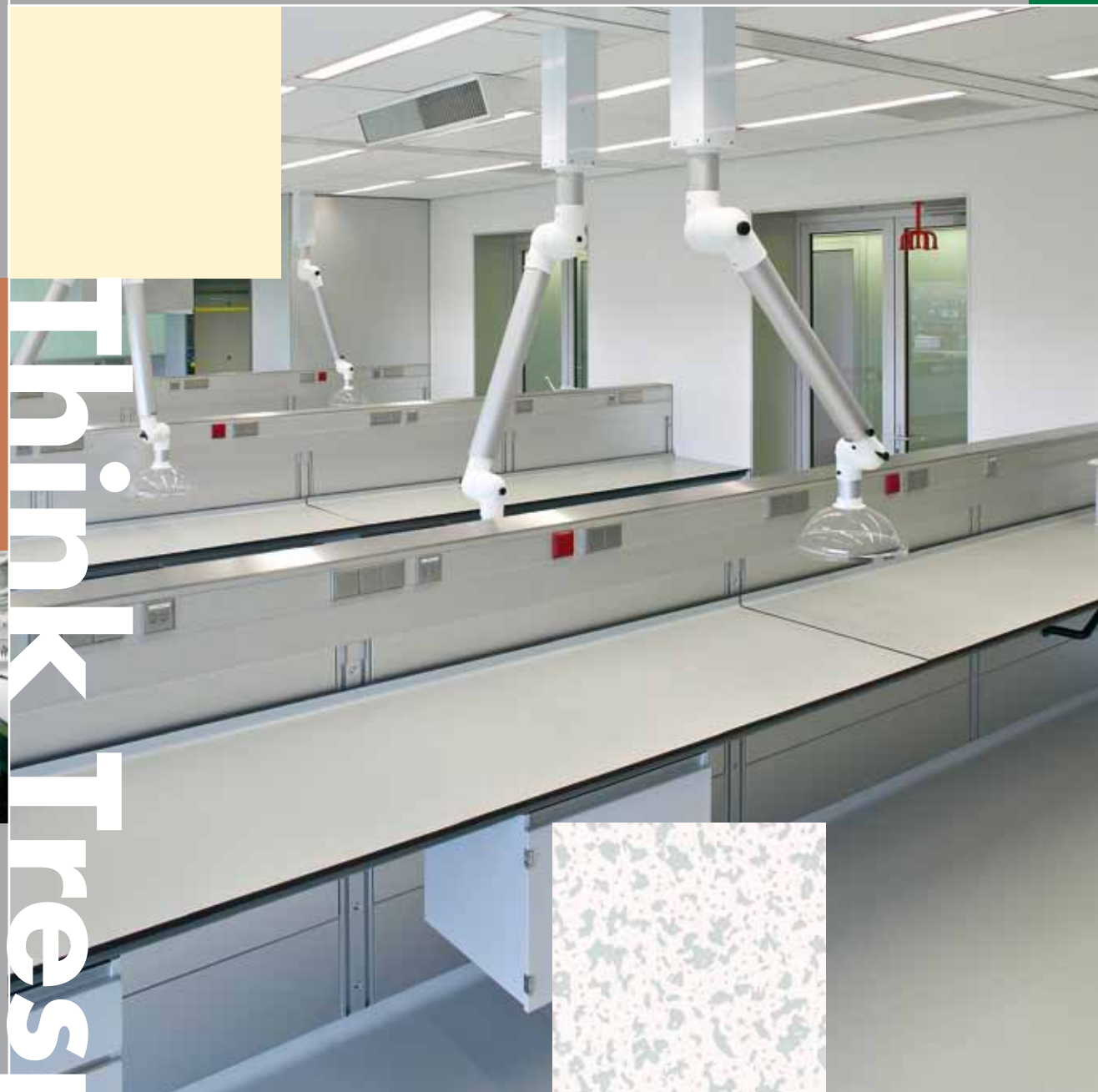

TRES|PA®

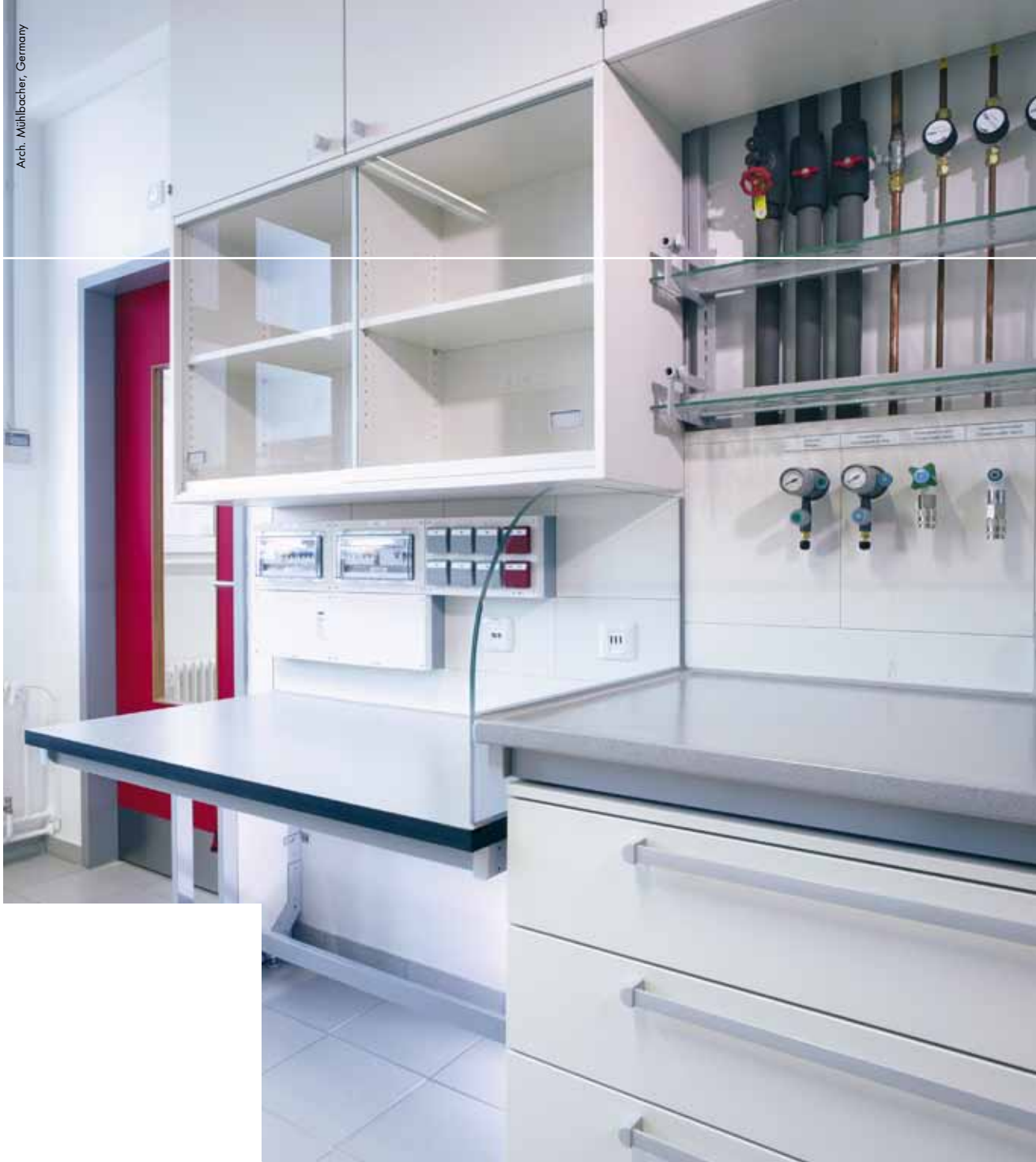

TRESPA TOPLAB^{PLUS} – FOR CLEAN AND STERILE SURFACES

Today's laboratory needs to meet the highest international standards. Therefore, selecting the right material for laboratory worktops and furniture is of crucial importance in an environment where tests and experiments must not be contaminated.

Hygiene, cleanliness, chemical and stain resistance are key requirements. Worktops in laboratories are especially vulnerable. The material used needs to be hygienic, easy to clean and maintain, durable and flexible.

Increasingly, there is also an additional requirement to find environmentally friendly solutions – to minimize ecological impact and maximize effectiveness.

Trespa TopLab^{PLUS} – the ideal and flexible solution

Trespa TopLab^{PLUS} panels meet the requirements of most laboratories while at the same time bringing a range of additional benefits. These benefits include chemical, wear and water resistance.

Trespa TopLab^{PLUS} is a self supporting flat panel, based on thermosetting resins, homogeneously reinforced with wood based fibers and manufactured under high pressure and at high temperatures. The panels have an integrated decorative surface using proprietary EBC technology. This process ensures that each panel is non-porous and resistant to a large number of aggressive chemicals, easy to clean, disinfect and maintain.

The surface is impermeable to most reagents used in all types of laboratories and is resistant to the effects of and wear and tear – making Trespa TopLab^{PLUS} ideal for use in labs utilized by different working groups, such as in educational institutions and industrial environments.

These properties make Trespa TopLab^{PLUS} highly suitable for use in medical laboratories and clean room conditions. It can be used in the most challenging environments for many years without losing its appearance or function.

Trespa TopLab^{PLUS} - a sustainable and green product

Environmental considerations play a significant role in the development and manufacture of Trespa TopLab^{PLUS}. For the production of Trespa TopLab^{PLUS} panels, Trespa uses a method for converting softwood from certified forests into an attractive, durable, moisture resistant and low maintenance material. Trespa TopLab^{PLUS} is an environmentally sound solution for laboratory worktop applications both today and tomorrow.

Trespa was one of the first manufacturers of construction products to map the totality of its products effects by means of a life-cycle analysis (LCA), right from the source of the raw materials to the end of the life-cycle.

Designing for the Future

- > Optimum panel sizes to reduce waste
- > Preferred panel thicknesses to reduce sub-frames
- > EBC surface technology to reduce cleaning and maintenance
- > Up to 70% wood based fibers made from sustainable resources
- > The majority of the raw materials used are made from sustainable resources
- > Ensuring emissions and discharge are within acceptable limits
- > Eliminating toxic and corrosive gases
- > Improving reduction of waste, storage and delay
- > Organizing efficient delivery and optimum routes






TRESPA TOPLAB^{PLUS} - A UNIQUE COMBINATION OF PROPERTIES

Antimicrobial performance

Thanks to its unique surface composition, Trespa TopLab^{PLUS} is non porous. Its anti-microbial properties are incorporated in the product without the use of coatings or additives. This means that these properties will remain active throughout the product's lifetime. Bacteria, molds and/or other microorganisms are unable to grow or penetrate the surface. Bacteria stains will dry relatively fast on the surface – and will not find any source of nutrition on the material. Independent tests by the British Industrial Microbiological Services Ltd (IMSL) show a reduction in bacteria of 99.99% after 24 hours.

Reduction of colony forming units (cm²)



on Trespa TopLab^{PLUS} surface

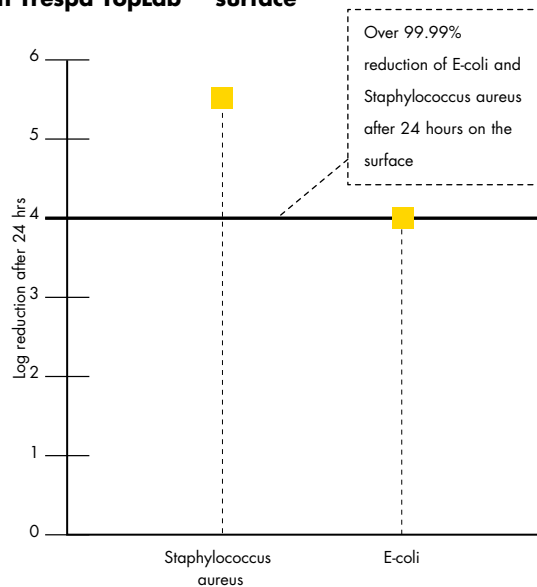

* Tested by using the Japanese Industrial Standard JIS Z 2801: 2000.

The photographs have been taken with the aid of a Scanning Electron Microscope. They clearly demonstrate the difference between a surface produced with Trespa EBC technology (Trespa TopLab^{PLUS}) and a surface of traditional melamine.

These microscopic pictures of the urethane-acrylic surface of Trespa TopLab^{PLUS} show that there is no difference between fresh material and material with which 10 years of aggressive cleaning is simulated.

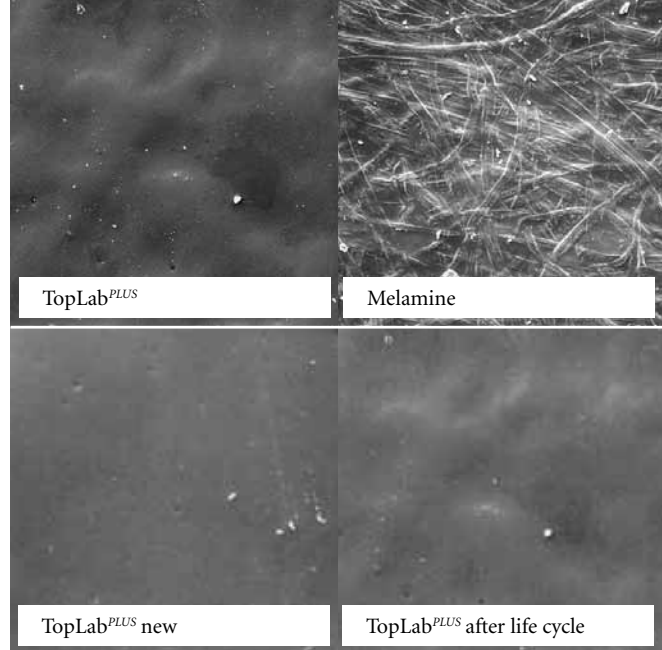

Arch. EGM, The Netherlands

Cleanability

In an environment where hygiene is key, Trespa TopLab^{PLUS} provides the best choice. Its surface is absolutely impervious to most materials used in biochemical and medical laboratories: radio-isotopes, human tissue and blood samples or bacteria. Biological or clinical test results are dependent on non-contamination. Trespa TopLab^{PLUS} provides a surface impermeable to most bacteria, molds or micro-organisms. Resistant to dyes and organic solvents, Trespa TopLab^{PLUS} is water-resistant and remains easy to clean or disinfect.

Thanks to the large size of Trespa TopLab^{PLUS} panels, the installation of large work tops with fewer joints is possible which fits the requirement for hygienic solutions. Moreover, for added hygiene, panels need not be glued onto a sub-structure.

Trespa TopLab^{PLUS} is self supporting from 1/2 in thickness and has a high load bearing ability.

Chemical resistance

Used extensively in chemical, analytical, micro-biological and educational laboratories world-wide, Trespa TopLab^{PLUS} is resistant to a large number of aggressive chemicals. Test results (see separate Trespa TopLab^{PLUS} test datasheet) show the panel's 24 hour resistance. Aggressive chemicals will not mark a Trespa TopLab^{PLUS} surface - if cleaned within this time period.

Durability

Trespa TopLab^{PLUS} is ideal for a multi-functional environment. Trespa TopLab^{PLUS} worktops are versatile and provide a tough and long-lasting surface that retains its appearance for many years. Used as part of a mobile and flexible environment, Trespa TopLab^{PLUS} delivers added strength to any laboratory or lecture room. The material's impact resistance makes it suitable for use in mobile furniture, i.e. trolleys.





TRESPA TOPLAB^{PLUS} – MEETING MARKET NEEDS

Design flexibility & fast transformation

Today's laboratory is a fast moving professional working environment. On the one hand it needs to house an array of testing and technical equipment as well as computers and their peripherals. On the other hand, there is a need for frequent changes to accommodate new tests or changes in operating requirements.

Trespa TopLab^{PLUS} facilitates the fast transformation of the laboratory. The material offers maximum design flexibility because it behaves in a similar way to hardwood. It can be machined and formed to meet the particular needs of the laboratory. Sinks, drainage holes, grooves and other accessories can be incorporated. Once installed, Trespa TopLab^{PLUS} can be easily adapted to accommodate changes in working practice. It can be recut and retrofitted with new taps, sinks or other equipment, without losing its exceptional performance characteristics and its appearance.

Pleasant working environment

A well appointed and attractive laboratory provides an efficient, effective and pleasant environment in which to work. Trespa TopLab^{PLUS} offers aesthetic benefits to the working environment. Panels are available in a range of colors to suit the specific needs of the facility's users.

Trespa TopLab^{PLUS} has a non-reflecting smooth surface. This makes it highly suitable as a multi-functional surface on which laboratory equipment, computers and general work, such as administration, analysis and research, can all be combined.

Trespa TopLab^{PLUS} is also available with an integrated decorative surface on both sides which widens the design possibilities and offers chemical and mechanical strength on those areas where it is needed on both sides, for example shelves and compartments. A workstation can be tailored to specific needs and still meet all the requirements of a laboratory.

Certification

Trespa TopLab^{PLUS} meets all the highest standards as laid down by leading national and international certifying authorities. Trespa TopLab^{PLUS} is GREENGUARD Indoor Air Quality Certified® and has achieved GREENGUARD Children & SchoolsSM Certification. Trespa TopLab^{PLUS} can also contribute to US Green Building Council LEED® credits. Trespa TopLab^{PLUS} holds the British Industrial Microbiological Services Ltd (IMSL) certification for anti-microbial performance, the certificate for hygiene from the German Institut für Krankenhaushygiene und Infektionskontrolle (IKI) and the German ISEGA certificate of conformity for contact with foodstuffs.





Sustainability

Trespa International strives to maximize the performance of its products, while minimizing their effects on the environment and on human health. Trespa was one of the first producers of panel materials to be certified according to ISO 14001, awarded by Lloyd's Register. The ISO 14001 standard describes the steps required for setting up, implementing, sustaining and improving a fully integrated environmental management system.


Trespa International supports its product-range with Life Cycle Analyses (LCAs) – studies which detail the environmental impact of each of their products throughout its life from “cradle to grave” – in terms of energy use, raw materials, etc.

Completely safe

The production of Trespa TopLab^{PLUS} is environmentally responsible and completely safe in use. Trespa TopLab^{PLUS} panels are based on thermosetting resin, homogeneously reinforced with up to 70% wood-based fibers. The majority of the used raw materials are made from sustainable resources.

A logical choice

Trespa TopLab^{PLUS} has been applied as the material of choice for many high standard laboratories around the world. With a strong emphasis on performance, cleanliness, environmental impact and aesthetics, Trespa TopLab^{PLUS} is the ideal surface solution for your laboratory.

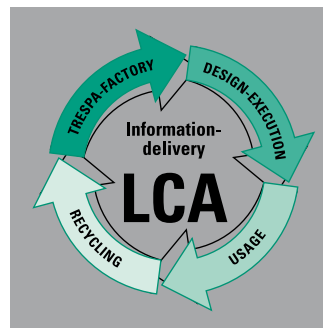





DISCLAIMER

General

These terms apply to the use of this document and such use automatically means that the other party agrees to these terms. The information provided by Trespa International B.V. ("Trespa") in this document is solely indicative. Trespa is unable to warrant the accuracy and completeness of this information. Trespa may change the information included in this document at any time and without further notice. Trespa's customers and third parties must ascertain that they have the most recent document (for the most recent version, please consult: www.trespa.com). No rights can be derived from the information provided; the use of the information is at the other party's risk and responsibility.

Trespa does not warrant that the information in this document is suitable for the purpose for which it is consulted by the other party. This document does not contain any design, structural calculation, estimate or other warranty or representation that customers and third parties may rely on. This document does not guarantee any properties of Trespa products. Colors used in Trespa's communications (including but not limited to printed matter) and in samples of Trespa's products may differ from the colors of the Trespa products to be supplied. Samples are not intended for use in product tests

and are not representative of characteristics of the Trespa products. Trespa's products and samples are produced within the specified color tolerances and the colors (of production batches) may differ, even if the same color is used. The viewing angle also influences the color perception. Metallics panels feature a surface whose color appears to change based on the direction from which it is viewed. The specified color stability and color specifications relate only to the decorative surface of the Trespa products, not to the core material and samples of the Trespa products. Trespa products are delivered ex-works with straight, sawn sides.

Customers and third parties must have a professional adviser inform them about (the suitability of) the Trespa products for all desired applications and about applicable laws and regulations. Trespa does not warrant the above.

The most recent version of the current delivery program and the Material Properties Datasheet can be found at www.trespa.info. Only the information in the most recent and valid Material Properties Datasheet should be used to select and provide advice regarding Trespa products. Trespa reserves the right to change (the specifications for) its products without prior notice.

Liability

Trespa is not liable (neither contractual nor non-contractual) for any damage arising from or related to the use of this document, except if and to the extent that such damage is the result of willful misconduct or gross negligence on the part of Trespa and/or its management. The limitation of liability applies to all parties affiliated with Trespa, including but not limited to its officers, directors, employees, affiliated enterprises, suppliers, distributors, agents, and representatives.

General conditions

All oral and written statements, offers, quotations, sales, supplies, deliveries and/or agreements and all related activities of Trespa are governed by the Trespa General Terms and Conditions of Sale (*Algemene verkooppvoorwaarden Trespa International B.V.*) filed with the Chamber of Commerce and Industry for Noord- en Midden-Limburg in Venlo (NL) on 11 April 2007 under number 24270677, which can be found on and downloaded from the Trespa website, www.trespa.com.

All oral and written statements, offers, quotations, sales, supplies, deliveries and/or agreements and all related work of Trespa North America, Ltd. are governed by the Trespa General Terms and Conditions of Sale, which can be found on and downloaded from the Trespa North

America Ltd. website, www.trespa.com/na. A copy of these general conditions of sale will be provided free of charge on request.

All general terms and conditions other than the conditions mentioned above are dismissed and do not apply, regardless of whether such terms and conditions are referred to on requests for offers, offer confirmations, stationery and/or other documents of the other party, even if Trespa does not expressly object to such terms and conditions.

Intellectual property

All intellectual property rights and other rights regarding the content of this document (including logos, text and photographs) are owned by Trespa and/or its licensors. Any use of the content of this document, including distribution, reproduction, disclosure, storage in an automated data file or the dispatch of such a file without Trespa's prior written consent is explicitly prohibited. © Trespa, Meteon, Athlon, TopLab, TopLab^{PLUS}, TopLab^{ECO-FIBRE}, Virtuon, Volkern, Trespa Essentials and Mystic Metallics are registered trademarks of Trespa.

Questions

Should you have any questions or comments, please do not hesitate to contact Trespa.

www.trespa.com

VISIT US

Trespa Design Centre New York

62 Greene Street (Ground Floor)
New York, NY 10012
United States of America
Tel.: 1-212-334-6888
Fax.: 1-866-298-3499
info.ny@trespa.com
www.trespa-ny.com

Trespa Design Centre Barcelona

Calle Ribera 5,
08003 Barcelona
Spain
Tel.: +34 (0) 93 295 4193
Fax.: +34 (0) 93 116 3300
info.barcelona@trespa.com
www.trespa-barcelona.com

CONTACT US

Trespa International B.V.

P.O. Box 110, 6000 AC Weert
Wetering 20, 6002 SM Weert
The Netherlands
www.trespa.com

Customer Service Desk EMEA Export

Tel.: +31 (0) 495 458 359 / 573
Fax: +31 (0) 495 458 383
infoexport@trespa.com

Trespa North America Ltd.

12267 Crosthwaite Circle
Poway, CA 92064
United States of America
Tel.: 1-800-487-3772
Fax: 1-858-679-9568
info.northamerica@trespa.com

Trespa China Co. Ltd.

Room 2604-05, HuaiHai Plaza
No. 1045 HuaiHai Road (central)
ShangHai 200031
P.R. China
Tel.: +86 (0) 21 6288 1299
Fax: +86 (0) 21 6288 1296
infochina@trespa.com

Customer Service Desk Asia/Pacific

Tel.: +86 (0) 21 5465 8388
Fax: +86 (0) 21 5465 6989
infoapac@trespa.com

Trespa UK Ltd.

Grosvenor House
Hollinswood Road
Central Park, Telford
Shropshire, TF2 9TW
United Kingdom
Tel.: +44 (0) 1952 290707
Fax: +44 (0) 1952 290101
info@trespa.co.uk




Your Trespa representative :

code V4651 version 3,0
date 08-18-2011 copies 0

Please visit www.trespa.com
for the most up to date version
of this document.

COLOR CARD TRESPA TOPLAB^{PLUS}

Uni colors

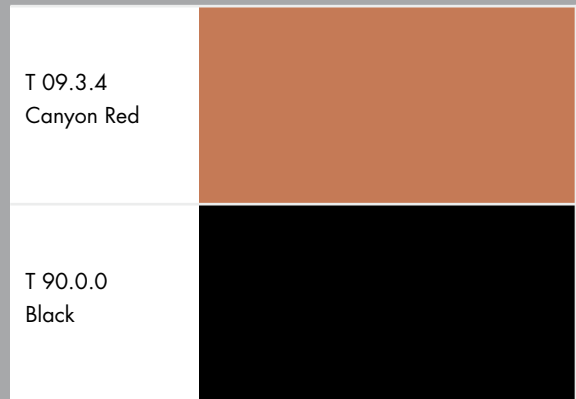

Speckles



The colors in this document are printed, and therefore, may vary slightly from the original Trespa panel colors with respect to gloss, color shades and surface texture. Original samples are available on request.

Disclaimer: All Trespa samples and all Trespa panels are produced within the specified tolerances. Samples and production panels do not originate from the same production batch. The color perception is by nature influenced by slight variations in pigments used as well as by the angle of observation.

TRESPA TOPLAB^{PLUS} - CHEMICAL RESISTANCE

(24 HOUR EXPOSURE)

Test procedure

The test was conducted by applying 5 drops of each reagent on the surface, covered with a watch glass (except those marked**). Chemicals marked ** were tested with a saturated cotton ball covered by a bottle. All chemicals were tested at room temperature for a period of 24 hours, rinsed off with water and evaluated.

Test results

No effect: No detectable stain, loss of gloss or change in work surface material.

Excellent: Slight stain or loss of gloss, but no change to the function, smoothness or life of the work surface material.

Good: A clearly discernible stain or loss of gloss, but no change to the function, smoothness or life of the work surface material.

Fair: Unacceptable staining or discernible deterioration or etching of the work surface material.

Failure: Severe stain or moderate deterioration, pitting cratering or etching of work surface material.

		No effect	Excellent	Good	Fair	Failure
Acids						
Acetic Acid	99%	■				
Acid Dichromate	5%	■				
Chromic Acid	60%	■				
Formic Acid	90%	■				
Hydrochloric Acid	10%	■				
Hydrochloric Acid	37%	■				
Hydrofluoric Acid	48%					■
Nitric Acid	20%	■				
Nitric Acid	30%		■			
Nitric Acid	65%			■		
Nitric Acid	70%			■		
Nitric Acid 65% : Hydrochloric Acid 37%	(1:3)	■				
Perchloric Acid	60%	■				
Phosphoric Acid	85%	■				
Sulphuric Acid	25%	■				
Sulphuric Acid	33%	■				
Sulphuric Acid	77%	■				
Sulphuric Acid	85%	■				
Sulphuric Acid	98%		■			
Sulfuric Acid 77% : Nitric Acid 70%	(1:1)			■		
Sulfuric Acid 85% : Nitric Acid 70%	(1:1)			■		
Bases						
Ammonium Hydroxide	28%	■				
Sodium Hydroxide	10%	■				
Sodium Hydroxide	20%	■				
Sodium Hydroxide	40%	■				
Sodium Hydroxide Flake		■				
Salts						
Copper Sulphate	10%	■				
Ferric(III)chloride	10%	■				
Potassium Iodite	10%	■				
Potassium Permanganate	10%	■				
Saturated Zinc Chloride		■				
Silver Nitrate	1%	■				
Sodium Chloride	10%	■				
Sodium Hypochlorite	13%	■				
Halogens						
Iodine (Crystals)			■			
Iodine Solution (0.1 N)			■			
Tincture of Iodine			■			
Organic Chemicals						
Cresol		■				
Dimethylformamide		■				
Formaldehyde	37%	■				
Furfural			■			

		No effect	Excellent	Good	Fair	Failure
Gasoline		■				
Hydrogen Peroxide	3%	■				
Phenol	90%	■				
Sodium Sulfide Saturated		■				
Solvents **						
Acetic Anhydride		■				
Acetone		■				
Acetonitrile		■				
Amyl Acetate		■				
Benzene		■				
Butyl Alcohol		■				
Carbon Tetrachloride		■				
Chloroform		■				
Dichlor Acetic Acid		■				
Dichloromethane		■				
Dioxane		■				
Diethyl Ether		■				
Ethylacetate		■				
Ethylalcohol		■				
Ethylene Glycol		■				
Methylalcohol		■				
Methylene Chloride		■				
Methylethylketone		■				
Methylisobutylketone		■				
Mono Chlorobenzene		■				
Napthelene		■				
n-Butyl Acetate		■				
Tetrahydrofurane		■				
n-Hexane		■				
Toluene		■				
Trichloroethylene		■				
Xylene		■				
Biological Stains						
Acridine Orange	1%	■				
Alizarin Complexone Dihydrate	1%	■				
Aniline Blue, water soluble	1%	■				
Basic Fuchsin	1%	■				
Carbol Fuchsin	1%	■				
Carmine	1%	■				
Congo Red	1%	■				
Gentian Violet (dye)	1%	■				
Eosin B	1%	■				
Giemsa Stain	1%	■				
Malachite Green Oxalate	1%	■				
Methyl Violet 2B	1%	■				
Methylene Blue	1%	■				
Safranin O	1%	■				
Sudan III	1%	■				
Wright Stain	1%	■				
Most conventional cleaning agents						
		■				



The chemicals in the above table include the 49 chemicals/concentrations set forth by SEFA 8 (Laboratory Casework) specifications as well as the main reagents from independent testing via Professional Service Industries/Pittsburgh Laboratory Division.

All information is based on our current state of knowledge. It is intended as information concerning our products and their application possibilities, and is therefore not intended as any form of guarantee with regard to any specific product characteristic. Test results differ per color.

Although the tests have been conducted according to the standard, it is recommended that users conduct their own tests: convince yourself that Trespa TopLab^{PLUS} is the only true multifunctional worktop!

Please check www.trespa.info for the latest version of the material properties and delivery program.